



PROmpter


Custom result reports for the Roche LightCycler PRO

welldone-prompter.github.io/PROmpter

Free & offline


Windows · one single file · no installation
no administrator rights · German & English

WORKFLOW From the instrument export to the finished report

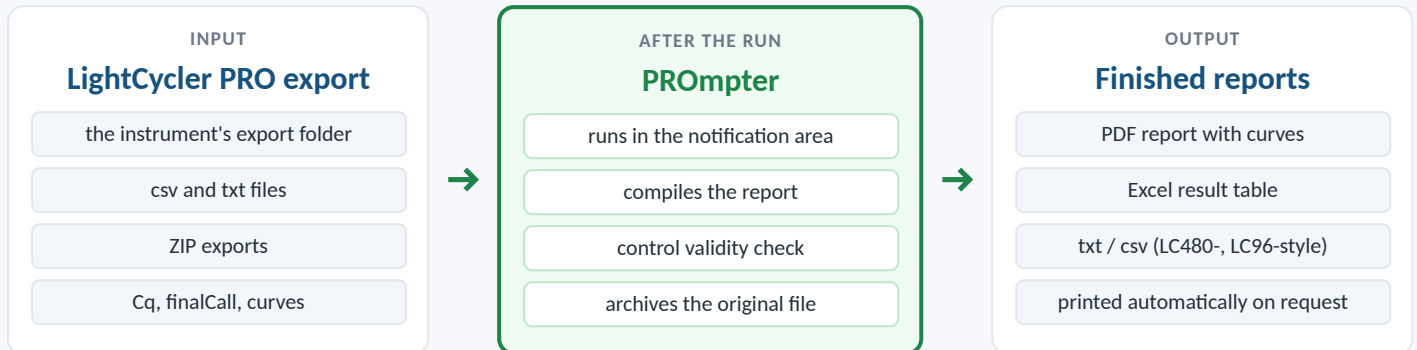

What PROmpter does

- **Automatic processing** – watches the LightCycler PRO export folder and runs in the background
- **PDF report** – findings, control validity check, plate overview and amplification curves
- **Excel table** – every result of a run, ready for further processing
- **txt/csv export** – in LC480 or LC96 format, or freely customisable
- **Customisable** – composition, naming, sorting and presentation of the reports
- **No recalculation** – Cq, finalCall and curves come unchanged from the LightCycler PRO
- **Bilingual** – interface and manual in German and English

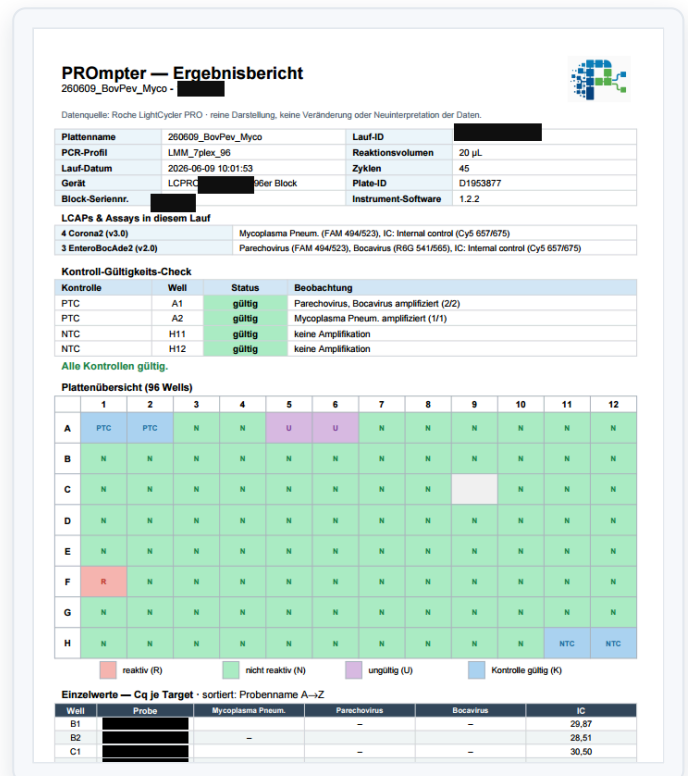

EXAMPLE OUTPUT

The finished PDF report – header data, control validity check, plate overview and the Cq values per target. Created automatically for every run (serial numbers redacted in this example).



Download PROmpter

welldone-prompter.github.io/PROmpter

Always the latest version – no installation, no administrator rights
Detailed manual (German / English) as a PDF on the same page



SISTER APP

WellDone!

Does the plate setup **before** the run:
pipetting scheme and LC PRO setup.

Sharing is explicitly encouraged – no sale

Feedback: welldone.pcr.setup@gmail.com

Private project, built with Claude Code

Important note: PROmpter is a helper tool, **not a medical device** and not validated. PROmpter is purely a presentation tool; the reports it creates are an interim evaluation, the professional release still happens manually. No liability is accepted for damage or incorrect results.

"Roche" and "LightCycler" are trademarks of their respective owners. PROmpter is an independent, private project, neither reviewed nor endorsed by Roche. Open-source libraries used: PySide6/Qt, ReportLab.